



Rebel Reliever

Patients diagnosed with mild to moderate medial or lateral OA, as well as patients who have OA in combination with ligament instability.

Features

- Rigid aircraft aluminium frame
- Patented adjustable Loadshifters
- Townsend Motion TM5+ Hinges
- Extension Stop Kit
- Exclusive Synergistic Suspension Strap
- Rotation control tibial shell contour
- Powdercoat finish



Large range of ordering options



Options

- Shell length
- Anti-migration strap padding
- Flexion stop kit
- A variety of colours

Ordering

- Standard OTS sizes S to XXL (tibia shell available in XS)
- "Customized" from 3 leg measurements
- Custom from a cast, digital leg scan, tracing, CAS measurement device, or Custom Configuration System (photos)


Townsend's Rebel Reliever OA knee brace has been clinically proven to reduce compressive forces in the knee by 36%. Other clinically proven benefits:

- Improved Walking Capacity
- Improved Muscle Function
- Reduced Compressive Forces
- Improved ADL Participation
- Reduced Symptoms
- Reduced Pain
- Improved Quality of Life
- Improved Balance Confidence

Size	6" above	M-L width	6" below
XS	12½ - 15½ inch	3 - 3½ inch	11 - 12½ inch
SM	15½ - 18½ inch	3½ - 4 inch	12¼ - 13¾ inch
MD	18½ - 21 inch	4 - 4½ inch	13¼ - 15 inch
LG	21 - 23½ inch	4½ - 5 inch	14¼ - 15¾ inch
XL	23½ - 25 inch	5 - 5½ inch	15 - 17 inch
XXL	25 - 28 inch	5½ - 6 inch	17 - 19 inch


Custom Fit From Three Measurements

Rebel Reliever braces ordered from 3 measurements accommodate patients with disproportional leg sizing. Standard braces made from only one leg measurement do not accommodate for patients who have disproportional legs.

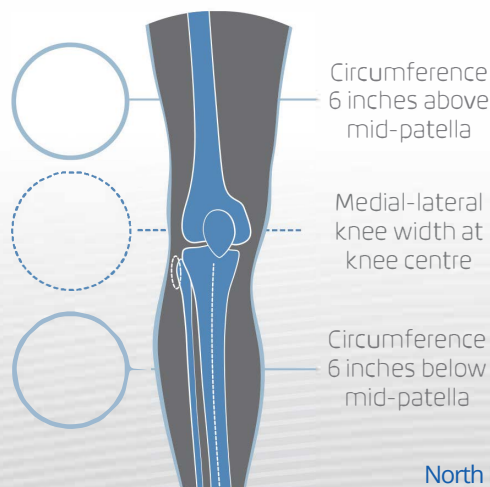

Corrective Force Technology

The uprights can be shortened and/or lengthened to shift the thigh shell angle up to 18 degrees. Increasing corrective force at the proximal end amplifies the affect of the three-point pressure system.

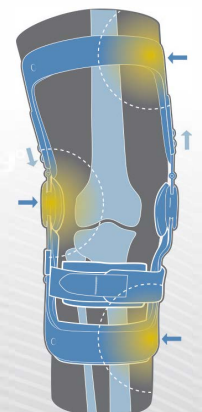

Thuasne UK Ltd

North Farm Road, Tunbridge Wells, TN2 3XF
Tel: 01295 257422 - Fax: 01295 257877 - Email: orders@thuasne.co.uk

Ordered by: _____ Phone # (_____) _____

BILLING: P.O. Number _____ Townsend Account # _____

Bill To: _____ Ship To: _____

Address: _____ Address: _____

City: _____ City: _____

County: _____ Postcode: _____ County: _____ Postcode: _____

Phone: (_____) _____ Fax: (_____) _____ Phone: (_____) _____ Fax: (_____) _____

Patient's Last Name: _____

Patient's First Name: _____

Male Female Age _____ Height _____ Weight _____

Compartment

Medial Compartment Lateral Compartment

Universal (select when ordering stock inventory)
Medial OA pads installed, Lateral OA pads included in kit

Leg Left Right

Shell Length Thigh shell extends 7" to 9" on all braces

Tibia 6 Inch 7 Inch 8 Inch

TM5+ Hinges Includes Extension Stops

Optional Flexion Stop Kit*

Standard Colors in stock

Argento Gray Textured Black

Non-stock Colors Brace may ship next business day

Bengal Silver Bengal Yellow Bengal White

Gloss Black Electric Blue Sky Blue Violet

Candy Green Sparkle Copper Sparkle Red

CS Package*

Size


X-Small Small Medium Large

X-Large XX-Large

Sizing Reference			
	6" above	M-L width	6" below
X-Small	NOT AVAILABLE	3" to 3.5"	11" to 12.5"
Small	15.5" to 18.5"	3.5" to 4"	12.25" to 13.75"
Medium	18.5" to 21"	4" to 4.5"	13.25" to 15"
Large	21" to 23.5"	4.5" to 5"	14.25" to 15.75"
X-Large	23.5" to 25"	5" to 5.5"	15" to 17"
XX-Large	25" to 28"	5.5" to 6"	17" to 19"

3 Measurements: "Customized" Assembly (No Added Charge)

If your patient has proportional leg sizing (see sizing reference, above) or if you are ordering a brace for stock inventory please select from the size options. However, if you would prefer to have Townsend customize the assembly of your patient's brace at no addition charge please provide leg measurements beside the illustration, below.

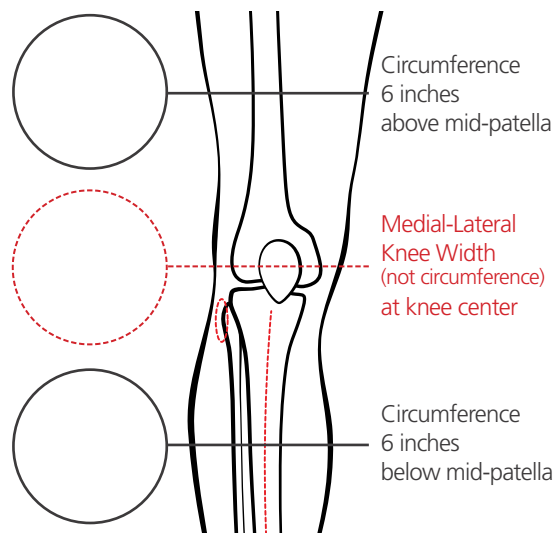

*Indicates additional charges apply

Thuasne UK Ltd

North Farm Road, Tunbridge Wells, TN2 3XF

Tel: 01295 257422 - Fax: 01295 257877

Email: orders@thuasne.co.uk