




THUASNE®

Available on
Prescription
NOW

BIFLEX Self Adjust

Adjustable compression wrap system for upper and lower limb with an average **cost saving of 34%** vs other compression wraps*

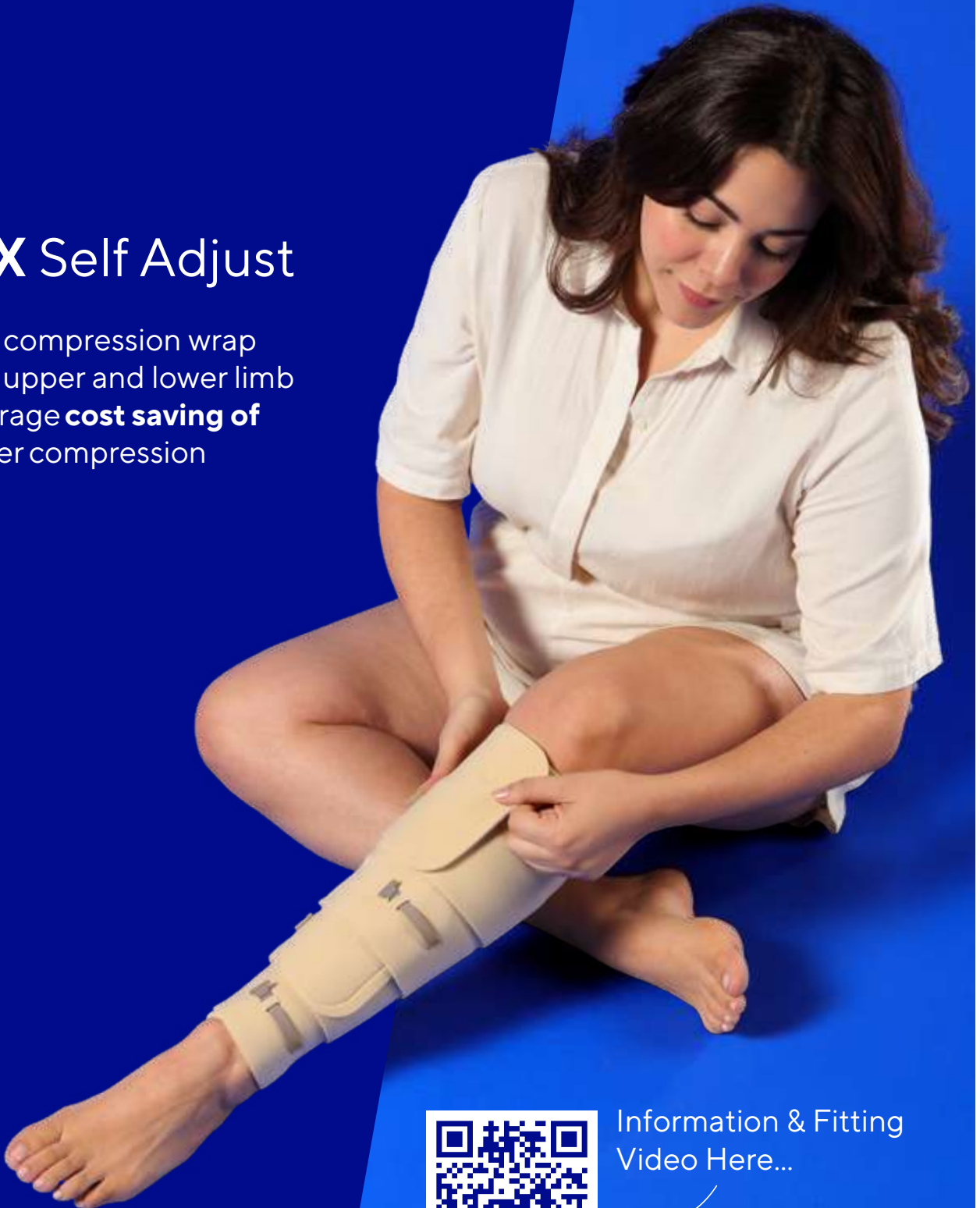

Information & Fitting
Video Here...




BIFLEX Self Adjust


Adjustable compression wrap system for upper and lower limb

Indicated for the management of oedema (venous, lymphatic and caused by injury), chronic venous insufficiency (including for the management of venous leg ulcers) and symptomatic treatment of lipoedema.


Patented, integrated pressure adjuster system on the calf and thigh garments. Delivering accurate 20-40mmHg pressures every time.*

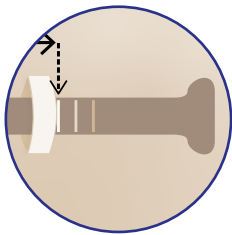


Comfortable and breathable open cell foam material, allowing moisture transfer from the skin.**

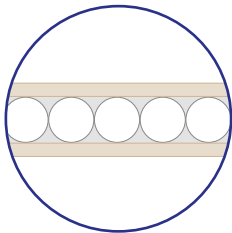

50% overlapping straps, mimicking compression bandaging and preventing oedema migration.

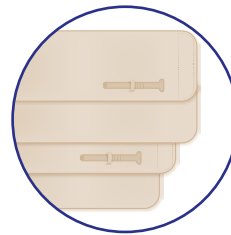

Easy donning and doffing - integrated liner (calf, thigh, arm), holds garment in position while applying the straps. Fewer straps for a simpler application.

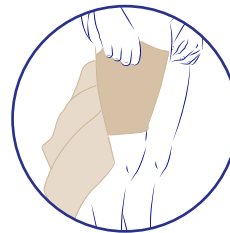


How to use the pressure adjuster ?

The adjuster comes with 3 levels of pressure: **L** low (light blue line), **M** medium (blue line) and **H** high (dark blue line).

Pressure is reached by pulling on the strap, which slides the adjuster until the outside of the loop is over the desired level of pressure.


If necessary, the segment pressure can be adjusted by adjusting the loops from bottom to top.

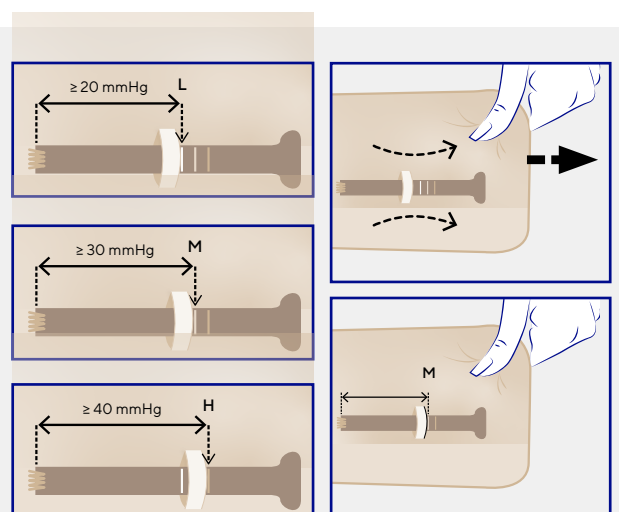

Size Guides


Calf

ℓ	Ankle-knee height
cB	Ankle circumference (narrowest point)
cC	Calf circumference (widest point)
cD	Circumference below the knee
N	Normal - $\ell < 36$ cm/ 14"
L	Long - $\ell > 36$ cm/ 14"

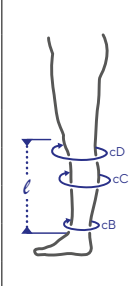
		cB	cC	cD
	1	16 - 24 cm 6.5 - 9.5"	24 - 35 cm 9.5 - 14"	25 - 36 cm 10 - 14"
	2	20 - 28 cm 8 - 11"	32 - 43 cm 12.5 - 17"	32 - 43 cm 12.5 - 17"
3	24 - 32 cm 9.5 - 12.5"	40 - 51 cm 15.5 - 20"	39 - 50 cm 15.5 - 19.5"	
	28 - 36 cm 11 - 14"	48 - 59 cm 19 - 23"	46 - 57 cm 18 - 22.5"	

Thigh

ℓ	Knee-upper thigh height
cE	Circumference above the knee
cG	Upper thigh circumference (widest point)
N	Normal - $\ell < 36$ cm/ 14"
L	Long - $\ell > 36$ cm/ 14"

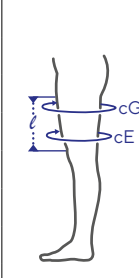
		cE	cG
	1	34 - 45 cm 13.5 - 18"	46 - 57 cm 18 - 22.5"
	2	41 - 53 cm 16 - 21"	53 - 65 cm 21 - 25.5"
3	49 - 60 cm 19.5 - 23.5"	61 - 72 cm 24 - 28.5"	
	56 - 68 cm 22 - 27"	68 - 80 cm 26.5 - 31.5"	

Knee

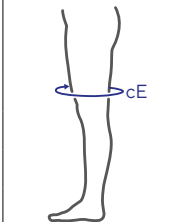		cE
	1	35 - 50 cm 14 - 19.5"
	2	50 - 65 cm 19.5 - 25.5"
3	65 - 85 cm 25.5 - 33.5"	

cE	Circumference above the knee
----	------------------------------

Foot

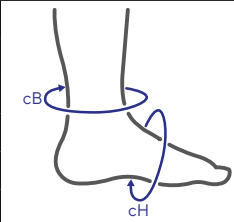		cB	cH
	1	< 36 cm < 14"	19 - 24 cm 7.5 - 9.5"
	2	< 36 cm < 14"	24 - 29 cm 9.5 - 11.5"
3	< 41 cm < 16"	29 - 34 cm 11.5 - 13.5"	


cB	Ankle circumference
cH	Instep circumference


Arm

cG	Axillary arm circumference
cD	Forearm circumference (widest point)
cC	Wrist circumference

Hand

cA	Hand circumference
----	--------------------

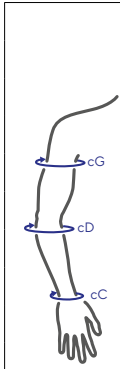		cG	cD	cC
	1	25 - 32 cm 10 - 12.5"	22 - 28 cm 9 - 11"	15 - 19 cm 6 - 7.5"
	2	29 - 36 cm 11.5 - 14"	24 - 30 cm 9.5 - 12"	16 - 20 cm 6.5 - 8"
3	33 - 40 cm 13 - 15.5"	26 - 32 cm 10 - 12.5"	17 - 21 cm 6.5 - 8"	
	37 - 44 cm 14.5 - 17.5"	28 - 34 cm 11 - 13.5"	17 - 21 cm 6.5 - 8"	
4	41 - 48 cm 16 - 19"	30 - 36 cm 12 - 14"	18 - 22 cm 7 - 8.5"	
	45 - 52 cm 18 - 20.5"	31 - 37 cm 12.5 - 14.5"	19 - 23 cm 7.5 - 9"	

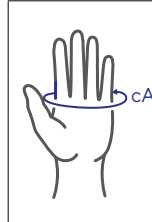		cA
	1	14 - 19 cm 5.5 - 7.5"
	2	19 - 24 cm 7.5 - 9.5"
3	24 - 29 cm 9.5 - 11.5"	

Scan QR Code for
Product & PIP Codes



WINGS
FOR YOUR HEALTH

BIFLEX Self Adjust

Adjustable compression wrap
system for upper and lower limb



Biflex Self Adjust is a medical device.
All the products mentioned on this document are class I medical devices, CE marked according to the European Regulation 2017/745 on medical devices. Please read carefully the instructions for use, indications and contraindications of the product.
Availability of these products might vary from a given country or region to another. Photos: Studio Caterin. Ref: 2306219 (06/2023).

 **THUASNE**
120 rue Marius AUFAN
92300 Levallois-Perret Cedex - France

Thuasne UK Ltd, Unit 4 Orchard Business Centre, North Farm Road, Tunbridge Wells, Kent, TN2 3XF.